

HEMOBLAST™ Bellows

HEMOSTATIC AGENT

Caution: Federal law restricts this device to sale on or by the order of a physician.




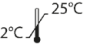









Note: A glossary of symbols used in the device labels is provided in Section 10.

1. Device Description	p.3	6. How Supplied	p.5
2. Indications for Use	p.3	7. Storage	p.5
3. Contraindications	p.3	8. Summary of Primary Clinical Studies	p.6
4. Warnings and Precautions	p.3	9. Potential Adverse Events	p.26
5. Directions for Use	p.5	10. Symbols Glossary	p.27

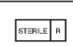

Manufacturer
Dilon Technologies France SAS
8, Allée Irène Joliot-Curie St-Priest
69800, France
www.dilon.com
Phone: +1 877 GO DILON
Fax: +1 757 269 4912

10. Symbols Glossary

Symbol	Standard	Reference #	Title	Definition
	ISO 15223-1	5.1.5	Batch code	Indicates the manufacturer's batch code so that the batch or lot can be identified.
	ISO 15223-1	5.1.6	Catalogue Number	Indicates the manufacturer's catalogue number so that the medical device can be identified.
	ISO 15223-1	5.1.4	Use-By date	Indicates the date after which the medical device is not to be used.
	ISO 15223-1	5.3.7	Temperature Limit Store between 2°C and 25°C (36-77°F)	Indicates the temperature limits to which the medical device can be safely exposed.
	ISO 15223-1	5.2.6	Do not re-sterilize	Indicates a medical device that is not to be re-sterilized.
	ISO 15223-1	5.2.8	Do not use if package is damaged and consult instructions for use	Indicates a medical device that should not be used if the package has been damaged or opened and that the user should consult the instructions for use for additional information.
	ISO 15223-1	5.4.2	Do not re-use	Indicates a medical device that is intended for one single use only.
	ISO 15223-1	5.4.3	Consult the instruction for use	Indicates the need for the user to consult the instructions for use.
	ISO 15223-1	5.2.4	Sterilized using irradiation	Indicates a medical device that has been sterilized using irradiation.
	ISO 15223-1	5.2.12	Double sterile barrier system	Indicates two sterile barrier systems.
	ISO 15223-1	5.6.3	Non-pyrogenic	Indicates a medical device that is non-pyrogenic.
	FDA Guidance "Alternative to Certain Prescription Device Labeling Requirements", issued 1/21/2000	N/A	Caution: Federal law restricts this device to sale on or by the order of a physician.	Requires prescription in the United States.
	ISO 15223-1	5.1.1	Manufacturer	Indicates the medical device Manufacturer

IFU-BQF02-US-V2
2023-02-20

Superiority in achieving hemostasis – 6 minutes

The secondary endpoint of superiority of HEMOBLAST™ in achieving hemostasis at 6 minutes was met in the overall population. As described in the Primary Endpoint section, the proportion of HEMOBLAST™ subjects achieving hemostasis at 6 minutes was 93.1% versus 73.5% for the G+T group (p=0.0001 for superiority).

Superiority in achieving hemostasis – 3 minutes

The secondary endpoint of superiority of HEMOBLAST™ in achieving hemostasis at 3 minutes was met in the overall population. As described in the Primary Endpoint section, the proportion of HEMOBLAST™ subjects achieving hemostasis at 3 minutes was 71.1% versus 45.8% for the G+T group (p=0.0001 for superiority).

In regard to surgical procedure and baseline TBS, the number of subjects in each surgical arm were found to be similar between the two treatment groups, as were the locations for the surgical procedure, the TBS tissue type, and the conventional procedures for hemostasis. The estimated size of the TBS was also found to be similar between treatment groups, with a mean of 5.4 cm² in the HEMOBLAST™ group and 5.8 cm² in the G+T group.

9. Potential Adverse Events

Risks possibly related to the use of hemostats similar to HEMOBLAST™ include:

- Adhesion formation;
- Allergy or anaphylaxis;
- Blockage of cardiopulmonary bypass system and cell saver devices;
- Compromised attachment of orthopedic implants;
- Creutzfeldt-Jakob disease (CJD) agent;
- Increased infection;
- Lack of efficacy;
- Nerve compression;
- Thrombosis or thromboembolism;
- Transmissible Spongiform Encephalopathies (TSE); and
- Viral disease transmission.

Hypersensitivity or allergic/anaphylactoid reactions may occur with HEMOBLAST™. Symptoms associated with allergic anaphylactic reactions include: flush, urticaria, pruritus, nausea, drop in blood pressure, tachycardia or bradycardia, dyspnea, severe hypotension, and anaphylactic shock.

These reactions may occur in patients receiving HEMOBLAST™ for the first time or may increase with repetitive applications of HEMOBLAST™. In the event of hypersensitivity reactions, discontinue administration of HEMOBLAST™. Mild reactions can be managed with antihistamines. Severe hypotensive reactions require immediate intervention using current principles of shock therapy.

1. Device Description

The HEMOBLAST™ Bellows consists of a powder composed of collagen, chondroitin sulfate, and thrombin. The powder is dry, sterilized, biocompatible, and non-pyrogenic. It is resorbed within 4 weeks. No preparation, mixing or heating is required.

HEMOBLAST™ Bellows is composed predominantly of highly purified porcine collagen with smaller amounts of bovine chondroitin sulfate and human pooled plasma derived thrombin. Each device contains a maximum of 1500 IU of thrombin. Plasma donations are from US plasma centers only. All individual donations of the plasma were tested for HBsAg, anti-HIV1/-HIV2, and anti-HCV and found to be negative. The plasma pools were tested and found to be non-reactive for HCV RNA, HBV DNA, and HIV1 RNA as determined by PCR (NAT). The product complies with the specifications of the manufacturer and World Health Organization (WHO).

The manufacturing procedures for the HEMOBLAST™ Bellows include processing steps designed to reduce the risk of viral transmission.

2. Indications for Use

HEMOBLAST™ Bellows is indicated in surgical procedures as an adjunct to hemostasis when control of minimal, mild, and moderate bleeding by conventional procedures is ineffective or impractical, except in neurosurgical, ophthalmic, and urological procedures.

3. Contraindications

- Do not inject HEMOBLAST™ Bellows into a vessel or tissue. There is a risk of allergic-anaphylactoid reaction and/or thromboembolic events which may be life-threatening.
- Do not apply HEMOBLAST™ Bellows in the absence of active blood flow, e.g., while the vessel is clamped or bypassed. Extensive intravascular clotting and even death may result.
- Do not use the HEMOBLAST™ Bellows for treatment of severe or extreme bleeding.
- Do not administer to patients with known allergies or hypersensitivity to materials of porcine or bovine origin.
- Do not use in the closure of skin incisions because it may interfere with the healing of the skin edges due to mechanical interposition of the powder.

4. Warnings and Precautions**Warnings:**

- The safety and effectiveness of HEMOBLAST™ Bellows has not been evaluated in ophthalmic, urologic, and neurologic procedures.
- Because HEMOBLAST™ Bellows is made from human blood, it may carry a risk of transmitting infectious agents, e.g., viruses, and theoretically, the Creutzfeldt-Jakob disease (CJD) agent. This also applies to any unknown or emerging viruses or other types of infections. The physician should discuss the risks and benefits of this product with the patient.
- HEMOBLAST™ Bellows contains chondroitin sulfate from bovine origin which is associated with a remote risk for Transmissible Spongiform Encephalopathies (TSE), which has been minimized in accordance with regulatory guidelines by a manufacturing process with demonstrated TSE inactivation capacity.
- The use of HEMOBLAST™ Bellows in laparoscopic or robotic assisted surgery procedures has not been studied and should not be used for acquiring hemostasis.
- HEMOBLAST™ Bellows is not intended as a substitute for meticulous surgical technique and the proper application of ligatures or other conventional procedures for hemostasis.
- When applied to a bleeding site, HEMOBLAST™ Bellows swells up to 60% within about 5 minutes. Caution should be used when applying HEMOBLAST™ Bellows in any cavities or closed spaces due to the swelling of the device.
- Do not attempt to trim the applicator tip.

- HEMOBLAST™ Bellows should not be used at the site of a valve replacement or repair as valvular dysfunction could occur.
- HEMOBLAST™ Bellows should not be applied at the site of a synthetic graft or patch implant due to potential for decreased effectiveness.
- The product should not be in contact with circulating cerebrospinal fluid (CSF).
- Do not implant in an infected or necrotic site or one that is likely to develop an infection.
- The safety and effectiveness of the combined use of the hemostatic agent with other devices has not been evaluated in controlled clinical trials.
- The material has not been tested on children and pregnant or lactating women.

Precautions:

- For single use only.
- Avoid direct contact between the nozzle or cannula tip and irrigation or bodily fluids during use to minimize the potential for clogging the tip during HEMOBLAST™ Bellows application.
- Do not re-sterilize. Consequences of re-sterilizing and/or re-using the hemostatic powder were not evaluated during the design. Therefore, the efficacy of the product may be reduced.
- Do not use HEMOBLAST™ Bellows in circumstances that result in negative peripheral venous pressure, e.g., patient positioning with head up may create negative pressure in venous structures above the diaphragm, as it may draw material into the vascular system potentially resulting in life-threatening thromboembolic events.
- Exposure to solutions containing alcohol, iodine, or heavy metals may cause the thrombin in hemostatic agent to be inactivated.
- Do not allow powder to enter into cell saver equipment, extracorporeal cardiopulmonary bypass circuits, or autologous blood salvage circuits.
- Do not use on bone surfaces where adhesives, such as methylmethacrylate or other acrylic adhesives, will be required to attach a prosthetic device.
- Avoid applying the HEMOBLAST™ Bellows device to open venous channels or sinuses of the liver parenchyma.
- HEMOBLAST™ Bellows should not be used for the primary treatment of coagulation disorders.
- Potential risk of thrombosis if absorbed systemically.
- Do not disrupt the HEMOBLAST™ Bellows clot complex by physical manipulation. This may cause clot disruption and subsequent bleeding.
- Do not use after the expiration date indicated on the package.
- Ensure that the packaging of each product is intact before use.
- Do not use if the package or the product is damaged.
- Do not use if the product has been exposed to extreme heat.
- The product must be manipulated by qualified personnel according to the general principles of sterility and pre-medication.
- The bellows should be handled with sterile powder-less gloves.
- There is a risk of human parvovirus B19 transmission. Human parvovirus B19 most seriously affects pregnant women or immune-compromised individuals. Symptoms of human parvovirus B19 include fever, drowsiness, chills, and runny nose followed about two weeks later by a rash and joint pain. Patients should report such symptoms to their physician if they appear.
- Long-term animal studies to evaluate the carcinogenic potential of HEMOBLAST™ Bellows or studies to determine the effect of HEMOBLAST™ Bellows on fertility have not been performed.
- It is not known whether HEMOBLAST™ Bellows can cause fetal harm when administered to pregnant women or can affect reproduction capacity.

The proportion of subjects in each treatment group achieving hemostasis at each time point was broken out by surgical arm, see Tables 19, 20 and 21. The HEMOBLAST™ group showed significantly higher rates of hemostasis at 3 and 6 minutes in the abdominal and orthopedic surgical arms. There was no significant difference between groups in the cardiothoracic arm, as shown in the tables below.

Table 19: Proportion Achieving Hemostasis for Cardiothoracic Surgical Arm

Time	All	HEMOBLAST™	G+T	P-value
3 minutes	31/58 (53.4%)	20/38 (52.6%)	11/20 (55.0%)	>0.9999
6 minutes	46/58 (79.3%)	31/38 (81.6%)	15/20 (75.0%)	0.7343
10 minutes	51/56 (91.1%)	34/37 (91.9%)	17/19 (89.5%)	>0.9999

Table 20: Proportion Achieving Hemostasis for Abdominal Surgical Arm

Time	All	HEMOBLAST™	G+T	P-value
3 minutes	46/90 (51.1%)	39/59 (66.1%)	7/31 (22.6%)	0.0001
6 minutes	74/90 (82.2%)	55/59 (93.2%)	19/31 (61.3%)	0.0003
10 minutes	80/88 (90.9%)	58/59 (98.3%)	22/29 (75.9%)	0.0015

Table 21: Proportion Achieving Hemostasis for Orthopedic Lower Extremity Surgical Arm

Time	All	HEMOBLAST™	G+T	P-value
3 minutes	74/90 (78.7%)	54/64 (87.1%)	20/32 (62.5%)	0.0082
6 minutes	89/90 (94.7%)	62/64 (100.0%)	27/32 (84.4%)	0.0037
10 minutes	92/90 (97.9%)	62/64 (100.0%)	30/32 (93.8%)	0.1135

Secondary endpoints

Secondary endpoints are reported for the overall population.

Mean preparation time

The preparation time for the HEMOBLAST™ group was found to be significantly shorter than the G+T group, with a mean of 0.37 minutes (22 seconds) for the HEMOBLAST™ group and a mean of 2.40 minutes (144 seconds) for the G+T group (-2.03 [-2.10, -1.86], p<0.0001). This shows that HEMOBLAST™ takes 2.03 minutes (122 seconds) less to prepare than the control agent.

Non-inferiority in achieving hemostasis – 3 min

The HEMOBLAST™ group met the secondary endpoint of non-inferiority in success achieving hemostasis at the 3 minute time point, with 71.1% of HEMOBLAST™ subjects achieving hemostasis versus 45.8% in the G+T group; a difference of 27.5% (14.0% to 40.9%, p<0.0001 for non-inferiority).

¹ HEMOBLAST™: muscle spasm, thoracic aorta dissection, hematoma, vocal cord paralysis, endoleak, pressure ulcer, decreased right ventricular function during cardiopulmonary bypass removal, arthralgia of right leg, recurrent shingles, gout flare, left thumb tenderness, insomnia, right shoulder acute bursitis, right bronchopleural fistula, blisters around tape covering incision; G+T: right atrium lead dislodgment, left arm weakness, femoral artery perforation.

Unanticipated Adverse Device Effects

There were no Unanticipated Adverse Device Effects.

The incidence of adverse events and serious adverse events were statistically the same in both experimental and control groups. The three main categories of potential adverse events identified in the pivotal study included thromboembolic events, wound healing complications primarily in the appearance of sternal dehiscence, and porcine collagen antibody titers.

2. Effectiveness Results

Primary endpoint

The primary efficacy endpoint of non-inferiority of HEMOBLAST™ relative to G+T for success at achieving hemostasis within 6 minutes was met, with 93.1% of the HEMOBLAST™ group achieving hemostasis at 6 minutes versus 73.5% of the G+T group, a difference of 19.5% (9.5% to 29.5%, p<0.0001 for non-inferiority) see Table 17. A 95% repeated confidence interval that accounts for the pre-specified stopping rule is calculated to be (7.1%, 31.9%), ruling out the non-inferiority margin of -10%.

Table 17: Primary Efficacy Endpoint with Full Efficacy Population

Time	HEMOBLAST™	G+T	Difference (95% CI)	Z-statistic	P-value
6 minutes	148/159 (93.1%)	61/83 (73.5%)	19.5% (9.5%, 29.5%)	5.7926	<0.0001

The confidence interval for the estimated difference in the probability of hemostasis at 6 minutes and the P value are computed using the Cochran-Mantel-Haenszel estimator stratified by surgical arms.

A summary of the proportion of each treatment group that achieved hemostasis at 3, 6, and 10 minutes is shown in Table 18. The HEMOBLAST™ group showed a higher proportion of patients achieving hemostasis at each time point assessed.

Table 18: Proportion of Each Treatment Group Achieving Hemostasis at 3, 6, and 10 Minutes

Time	All	HEMOBLAST	G+T
3 minutes	151/242 (62.4%)	113/159 (71.1%)	38/83 (45.8%)
6 minutes	209/242 (86.4%)	148/159 (93.1%)	61/83 (73.5%)
10 minutes	223/238 (93.7%)	154/158 (97.5%)	69/80 (86.3%)

Numbers are n/N (percent).

The study was stopped early for efficacy, per the IDMC recommendation, based on the pre-specified stopping rules.

5. Directions for Use

Inspect the integrity of the product packaging prior to use. If damaged, do not use. The instructions must be read before use.

1. Peel open the trays/blisters to remove contents with attention to sterile procedures.
2. Remove cap using a twisting motion.
3. Attach the 10 cm nozzle extension to the nozzle when additional nozzle length is needed.
4. Blot excess blood from the target bleeding site with gauze/pad or suction so the hemostatic implant material may be applied directly to the source of bleeding. The wound surface should be as dry as possible just before application.
5. Apply hemostatic implant material to the source of bleeding by compressing the bellows. Enough implant material should be applied to cover the entire source of bleeding.
6. Immediately use a wet laparotomy pad to hold the hemostatic material at the target bleeding site against the bleeding surface using wound appropriate pressure to conform the powder to the source of bleeding.
7. Maintain the hemostatic material at the target bleeding site for approximately three minutes.
8. Gently lift the laparotomy pad and inspect the area.
9. If hemostasis has not been achieved, repeat steps 4-7 or use an alternate method of hemostasis treatment.
10. Discard any unused product after opening.

Notes:

- Multiple layers of the hemostatic agent may be applied onto the wound to cover the bleeding surface.
- When applied with the nozzle tip 2-3 cm from the target bleeding site, HEMOBLAST™ Bellows can treat a surface area up to 50 cm².
- Do not disrupt the HEMOBLAST™ Bellows clot complex by physical manipulation.
- Once bleeding has ceased, excess HEMOBLAST™ Bellows powder not incorporated in the hemostatic clot should be removed by gentle irrigation.
- Once hemostasis is achieved, excess HEMOBLAST™ Bellows should be removed from the site of application by irrigation and aspiration particularly when used in and around the posterior fossa, foramina of bone, areas of bony confine, the spinal cord, and/or the optic nerve and chiasm.
- HEMOBLAST™ Bellows should not be used for controlling post-partum bleeding or menorrhagia.
- In urological procedures, HEMOBLAST™ Bellows should not be left in the renal pelvis or ureters to eliminate the potential foci for calculus formation.

6. How Supplied

HEMOBLAST™ Bellows is supplied as a bellows pre-loaded with hemostatic powder. The bellows contains 1.65g of implant material. A 10 cm nozzle extension is included to assist with the application of the powder to active bleeding sites where a slightly longer tip is desired by the surgeon. The product is sterilized using gamma-sterilization and provided in double-packaging.

7. Storage

- The product is sterile unless the package has been opened, damaged, or otherwise contaminated.
- The product must be stored between 2°C and 25°C (36-77°F) in a dry environment.
- The product can be used throughout the duration of the operation without any known or reported hazards.
- The product must be discarded after use.

8. Summary of Primary Clinical Studies

Biom'up performed clinical studies to establish a reasonable assurance of safety and effectiveness of general surgical procedures with Hemoblast Bellows as an adjunct to hemostasis when control of minimal, mild, and moderate bleeding by conventional procedures is ineffective or impractical, except in neurosurgical, ophthalmic, and urological procedures in the US under IDE # G150037 and # G160063. Data from these clinical studies were the basis for the PMA approval decision. A summary of the clinical studies is presented below.

Pilot Clinical Study:

A. Study Design

This was a prospective, multicenter, single-arm, pilot, clinical investigation to collect data to support the use of the Surface Bleeding Severity Scale ("SBSS"), and to collect initial safety and efficacy data on HEMOBLAST™ Bellows in clinical use. The pilot study enrolled 31 subjects; 27 subjects were included in the safety population and 24 subjects were included in the effectiveness population. Subjects enrolled into the study were required to have a target bleeding site with SBSS score of 1 (minimal bleeding), 2 (mild bleeding), or 3 (moderate bleeding). This 6-week study was a single arm study which included patients undergoing only orthopedic and abdominal surgeries with associated bleeding sites; the was no cardiothoracic arm.

The bleeding scale elements used in the SBSS are provided in the table below.

Table 1: Surface Bleeding Severity Scale

Surface Bleeding Severity Scale Score	0	1	2	3	4	5
Verbal Descriptor	None	Minimal	Mild	Moderate	Severe; not immediately life-threatening	Extreme; immediately life-threatening
Visual Descriptor	Dry	Oozing	Pooling	Flowing	Streaming	Gushing
Expected Intervention(s)	None	Manual pressure, cautery, adjuvant hemostat(s)	Manual pressure, cautery, suture, adjuvant hemostat(s)	Manual pressure, cautery, suture, adjuvant hemostat(s)	Manual pressure, cautery, suture, staples, tissue repair	Manual pressure, cautery, suture, staples, tissue repair
Maximum Expected ACS-ATLS Shock Risk Class	1	1	1	2	3	4
<p>ACS-ATLS Shock Risk Class: 1 – involves up to 15% of blood volume; typically no change in vital signs and fluid resuscitation is necessary. Class 2 – involves 15-30% of total blood volume; patient is often tachycardic with a narrowing of the difference between the systolic and diastolic blood pressures; the body attempts to compensate with peripheral vasoconstriction; skin may start to look pale and be cool to the touch; volume resuscitation with crystalloids is all that is typically required; blood transfusion is not typically required. Class 3 – involved loss of 30-40% of circulating blood volume; patient's blood pressure drops; heart rate increases, peripheral hypoperfusion worsens; fluid resuscitation with crystalloid and blood transfusion are usually necessary. Class 4 – involves loss of > 40% of circulating blood volume; the limit of the body's compensation is reached and aggressive resuscitation is required to prevent death.</p> <p>References: 1. Spotnitz WD, Zielske D, Centis V et al. The SPOT GRADE: a new method for reproducibly quantifying surgical wound bleeding. Spine (Phila Pa 1976). 2017 Oct 10. doi: 10.1097/BRS.0000000000002447. [Epub ahead of print]. 2. Kortbeek JB, Al Turki SA, Ali J et al. Advanced trauma life support. 8th edition, the evidence for change. J Trauma 2008;64:1638-50.</p>						

Adverse Event Type*	All	HEMOBLAST™	G+T
Hematoma	3/258 (1.2%)	2/175 (1.1%)	1/83 (1.2%)
Infection	4/258 (1.6%)	3/175 (1.7%)	1/83 (1.2%)
Non-healing	2/258 (0.8%)	1/175 (0.6%)	1/83 (1.2%)
Seroma	6/258 (2.3%)	6/175 (3.4%)	0/83 (0.0%)
Other ⁴	1/258 (0.4%)	0/175 (0.0%)	1/83 (1.2%)
Other ⁵	11/258 (4.3%)	9/175 (5.1%)	2/83 (2.4%)

* Open text field, like responses were pooled.

¹HEMOBLAST™: Pancolitis, diarrhea positive for clostridium difficile, dental abscess; G+T: tricuspid valve vegetation

²HEMOBLAST™: right leg pain, back pain, right ear pain, left hip pain, right shoulder; G+T: sore shoulder, pain at epidural site and back pain, abdominal pain

³Intraoperative TEE demonstrated internal jugular vein clot

⁴Wound irritation from one stitch

⁵HEMOBLAST™: muscle spasm, thoracic aorta dissection, hematoma, vocal cord paralysis, endoleak, pressure ulcer, decreased right ventricular function during cardiopulmonary bypass removal, arthralgia of right leg, recurrent shingles, gout flare, left thumb tenderness, insomnia, right shoulder acute bursitis, right bronchopleural fistula, blisters around tape covering incision; G+T: right atrium lead dislodgment, left arm weakness, femoral artery perforation

Non-serious adverse events

Table 16 lists non-serious AEs occurring in 5% or more of all patients, or in one of the treatment arms.

Table 16: Number of Subjects Experiencing Each Type of Non-serious Adverse Event

Adverse Event Type*	All	HEMOBLAST	G+T
Abnormal Bloodwork	18/258 (7.0%)	13/175 (7.4%)	5/83 (6.0%)
Anemia	19/258 (7.4%)	10/175 (5.7%)	9/83 (10.8%)
Arrhythmia	23/258 (8.9%)	16/175 (9.1%)	7/83 (8.3%)
Constipation/Ileus	17/258 (6.6%)	12/175 (6.9%)	5/83 (6.0%)
Fluid Overload	13/258 (5.0%)	8/175 (4.6%)	5/83 (6.0%)
Infection (Non wound Related)	13/258 (5.0%)	7/175 (4.0%)	6/83 (7.2%)
Nausea	18/258 (7.0%)	15/175 (8.6%)	3/83 (3.6%)
Pain	36/258 (14.0%)	25/175 (14.3%)	11/83 (13.3%)
Wound Related	16/258 (6.2%)	12/175 (6.9%)	4/83 (4.8%)
Other ¹	11/258 (4.3%)	9/175 (5.1%)	2/83 (2.4%)

* Open text field, like responses were pooled.

Adverse Event Type*	All	HEMOBLAST™	G+T
<i>Yeast</i>	2/258 (0.8%)	1/175 (0.6%)	1/83 (1.2%)
<i>Other¹</i>	4/258 (1.6%)	3/175 (1.7%)	1/83 (1.2%)
Lactic Acidosis	1/258 (0.4%)	1/175 (0.6%)	0/83 (0.0%)
Lethargy	1/258 (0.4%)	0/175 (0.0%)	1/83 (1.2%)
Leukocytosis	3/258 (1.2%)	3/175 (1.7%)	0/83 (0.0%)
Nausea	18/258 (7.0%)	15/175 (8.6%)	3/83 (3.6%)
Orthostasis	1/258 (0.4%)	0/175 (0.0%)	1/83 (1.2%)
Pain	36/258 (14.0%)	25/175 (14.3%)	11/83 (13.3%)
<i>Mouth/Throat</i>	8/258 (3.1%)	5/175 (2.9%)	3/83 (3.6%)
<i>Post-Operative</i>	28/258 (10.9%)	21/175 (12.0%)	7/83 (8.4%)
<i>Other²</i>	8/258 (3.1%)	5/175 (2.9%)	3/83 (3.6%)
Paresthesia	1/258 (0.4%)	1/175 (0.6%)	0/83 (0.0%)
Pericarditis	1/258 (0.4%)	1/175 (0.6%)	0/83 (0.0%)
Pneumothorax	2/258 (0.8%)	0/175 (0.0%)	2/83 (2.4%)
Pruritis	4/258 (1.6%)	3/175 (1.7%)	1/83 (1.2%)
Respiratory Insufficiency	9/258 (3.5%)	6/175 (3.4%)	3/83 (3.6%)
TBS Re-bleed	5/258 (1.9%)	1/175 (0.6%)	4/83 (4.8%)
Thrombocytopenia	2/258 (0.8%)	2/175 (1.1%)	0/83 (0.0%)
Thromboembolic Event	4/258 (1.6%)	4/175 (2.3%)	0/83 (0.0%)
<i>Thrombosed AV Fistula</i>	1/258 (0.4%)	1/175 (0.6%)	0/83 (0.0%)
<i>Deep Vein Thrombosis</i>	2/258 (0.8%)	2/175 (1.1%)	0/83 (0.0%)
<i>Other³</i>	1/258 (0.4%)	1/175 (0.6%)	0/83 (0.0%)
Urinary Retention/Oliguria	7/258 (2.7%)	4/175 (2.3%)	3/83 (3.6%)
Vaginal Discharge	1/258 (0.4%)	0/175 (0.0%)	1/83 (1.2%)
Vomiting	5/258 (1.9%)	4/175 (2.3%)	1/83 (1.2%)
Wound Related	16/258 (6.2%)	12/175 (6.9%)	4/83 (4.8%)
<i>Dehiscence</i>	1/258 (0.4%)	1/175 (0.6%)	0/83 (0.0%)

1. Clinical Inclusion and Exclusion Criteria

A subject must have met all the following preoperative inclusion criteria to be enrolled into the investigation:

- Subject is undergoing elective open abdominal or orthopedic lower extremity surgery;
- Subject or an authorized legal representative is willing and able to give prior written informed consent for investigation participation;
- Subject on antiplatelets, including aspirin, will discontinue medication at least 10 days prior to surgery; and
- Subject is 21 years of age or older.

A subject had to meet all of the following intraoperative inclusion criteria to be enrolled into the investigation:

- Subject does not have an active or suspected infection at the surgical site;
- Subject in whom the Investigator is able to identify a target bleeding site (TBS) for which any applicable or conventional means for hemostasis are ineffective or impractical; and
- Subject has a TBS with SBSS score of 1, 2, or 3.

A subject had to meet all of the following preoperative exclusion criteria to be enrolled into the investigation:

- Subject is undergoing a laparoscopic, thoracoscopic, or robotic surgical procedure;
- Subject is undergoing a spinal surgical procedure;
- Subject is undergoing a neurologic surgical procedure;
- Subject is undergoing an emergency surgical procedure;
- Subject is pregnant, planning on becoming pregnant during the follow-up period, or actively breast-feeding;
- Subject has a clinically significant coagulation disorder or disease, defined as a platelet count <100,000 per microliter and/or International Normalized Ratio > 1.5 within 4 weeks of surgery;
- Subject had chronic corticosteroid use within 2 weeks prior to surgery;
- Subject receiving intravenous heparin or oral Coumadin within 24 hours of surgery;
- Subject has an active or suspected infection at the surgical site;
- Subject has had or has planned any organ transplantation;
- Subject has a known sensitivity or allergy to bovine and/or porcine substance(s) or any other component(s) of the hemostatic agent;
- Subject has ASA classification of > 4;
- Subject has a life expectancy of less than 3 months;
- Subject has a known psychiatric disorder, which in the opinion of the Principal Investigator, would preclude the subject from completing this clinical study;
- Subject has documented severe congenital or acquired immunodeficiency;
- Subject has religious or other objections to porcine or bovine components;
- Subject in whom the investigational device will be used at the site of a cemented or uncemented porous coated joint implant;
- Subject is currently participating or has participated in another clinical trial within the past 30 days and is receiving/has received an investigational drug, device, or biologic agent; and
- Subject is not appropriate for inclusion in the clinical trial, per the medical opinion of the Principal Investigator.

2. Follow-up Schedule

All enrolled subjects were evaluated in the immediate postoperative period (within 24 hours after surgery) and at 6 weeks postoperatively (visit 4-8 weeks postoperatively).

3. Clinical Endpoints

Primary endpoint:

The primary endpoint of this clinical investigation is the mean paired Kappa statistic for the assignment of SBSS scores by 2 Investigators.

Secondary endpoints:

Secondary endpoints of this clinical investigation consist of:

- Proportion of subjects achieving hemostasis within 6 minutes of HEMOBLAST™ Bellows application;
- Proportion of subjects achieving hemostasis within 10 minutes of HEMOBLAST™ Bellows application;
- Proportion of subjects achieving hemostasis within 3 minutes of HEMOBLAST™ Bellows application; and
- Incidence of adverse events through final follow-up.

B. Accountability

A total of 31 subjects met preoperative eligibility criteria, including roll-in subjects. Four subjects were excluded from the study intraoperatively due to not meeting intraoperative inclusion criteria, resulting in a total of 27 enrolled subjects, including the roll-in subjects. These 27 subjects are considered the safety analysis population.

There was one roll-in subject at each site and roll-in subjects are not included in the efficacy analysis. Therefore, 24 subjects enrolled in the study are considered the efficacy analysis population. The exploratory analysis population includes 25 subjects total (the efficacy analysis population plus one subject who met pre-operative eligibility requirements but failed intra-operative eligibility criteria). Enrollment was fairly balanced across all 3 sites, with 9 subjects being enrolled at Site 2 and 8 and 10 subjects enrolled at Site 6 and Site 8, respectively. Site 2 had 3 subjects that were excluded intraoperatively, Site 8 had 1 subject that was excluded intraoperatively, while no subjects were excluded intraoperatively at Site 6.

Table 2: Number of Subjects per Site

Population	All	Site 2	Site 6	Site 8
Total: Preoperative eligible + roll-in subjects	31	12	8	11
Intraoperative ineligible subjects	4	3	0	1
SAFETY ANALYSIS POPULATION: Preoperative AND Intraoperative eligible (ENROLLED and includes roll-in subjects)	27	9	8	10
Roll-in subjects	3	1	1	1
EFFICACY ANALYSIS POPULATION: Preoperative AND Intraoperative eligible (ENROLLED and EXCLUDES roll-in subjects)	24	8	7	9

Adverse Event Type*	All	HEMOBLAST™	G+T
Anemia	19/258 (7.4%)	10/175 (5.7%)	9/83 (10.8%)
Anxiety	1/258 (0.4%)	0/175 (0.0%)	1/83 (1.2%)
Arrhythmia	23/258 (8.9%)	16/175 (9.1%)	7/83 (8.3%)
Atelectasis	7/258 (2.7%)	4/175 (2.3%)	3/83 (3.6%)
Bile Leak	1/258 (0.4%)	1/175 (0.6%)	0/83 (0.0%)
Constipation/Ileus	17/258 (6.6%)	12/175 (6.9%)	5/83 (6.0%)
Dehydration	1/258 (0.4%)	0/175 (0.0%)	1/83 (1.2%)
Delirium	2/258 (0.8%)	1/175 (0.6%)	1/83 (1.2%)
Depression	1/258 (0.4%)	1/175 (0.6%)	0/83 (0.0%)
Dislocation	2/258 (0.8%)	1/175 (0.6%)	1/83 (1.2%)
Diarrhea	3/258 (1.2%)	3/175 (1.7%)	0/83 (0.0%)
Dizziness	3/258 (1.2%)	2/175 (1.1%)	1/83 (1.2%)
Epistaxis	1/258 (0.4%)	1/175 (0.6%)	0/83 (0.0%)
Fever	4/258 (1.6%)	2/175 (1.1%)	2/83 (2.4%)
Fluid Overload	13/258 (5.0%)	8/175 (4.6%)	5/83 (6.0%)
<i>Hypervolemia</i>	1/258 (0.4%)	0/175 (0.0%)	1/83 (1.2%)
<i>Edema</i>	6/258 (2.3%)	3/175 (1.7%)	3/83 (3.6%)
<i>Pleural Effusion</i>	6/258 (2.3%)	5/175 (2.9%)	1/83 (1.2%)
Fracture	1/258 (0.4%)	1/175 (0.6%)	0/83 (0.0%)
Hypertension	3/258 (1.2%)	2/175 (1.1%)	1/83 (1.2%)
Hypotension	3/258 (1.2%)	2/175 (1.1%)	1/83 (1.2%)
Infection (Non-wound Related)	13/258 (5.0%)	7/175 (4.0%)	6/83 (7.2%)
<i>Bacteremia</i>	3/258 (1.2%)	2/175 (1.2%)	1/83 (1.2%)
<i>Drain</i>	1/258 (0.4%)	1/175 (0.6%)	0/83 (0.0%)
<i>Phlebitis</i>	1/258 (0.4%)	0/175 (0.6%)	1/83 (1.2%)
<i>Respiratory</i>	3/258 (1.2%)	1/175 (0.6%)	2/83 (2.4%)
<i>Urinary Tract</i>	4/258 (1.6%)	1/175 (0.6%)	3/83 (3.6%)

D. Safety and Effectiveness Results

1. Safety results

The analysis of safety was based on the pivotal study cohort of 175 patients available for the 6 week evaluation. The key safety outcomes and adverse effects for this study are presented below.

In summary, at the time of the interim analysis, the HEMOBLAST™ group had a lower rate of subjects experiencing adverse events, although this did not reach statistical significance (48.0% vs. 56.6%, p=0.1516).

Serious adverse events

The proportion of patients experiencing a serious adverse event was comparable between treatment groups (10.9% for HEMOBLAST™ vs. 13.3% for G+T, p=0.5415). The two groups were found to have a similar number of subjects experiencing an adverse event related to the device.

There were 3 possible device-related adverse events reported in 3 subjects receiving HEMOBLAST™:

- Subject reported symptoms indicating possible systemic inflammatory response syndrome. This event was classified as possibly related to HEMOBLAST™; the event resolved with no sequelae.
- Subject experienced a generalized skin reaction with hives 11 days postoperatively; the event resolved with no sequelae. This was deemed as possibly related to HEMOBLAST™ and the CEC adjudicated this event as a Serious Adverse Device Effect (SADE).
- Subject was diagnosed with a pulmonary embolus 10 days postoperative; the event resolved with no sequelae. This was deemed as possibly related to HEMOBLAST™, but probably related to the investigational procedure.

There were 4 possible device-related adverse events reported in 4 subjects receiving G+T:

- Subject experienced acute blood loss anemia 4 days postoperative and was transfused with packed RBCs; the event resolved with no sequelae. This event was classified as possibly related to G+T.
- Subject developed a period of hypotension following surgery; it resolved with no sequelae after an infusion of albumin. This event was classified as possibly related to G+T.
- Subject had re-bleeding of the target bleeding site treated with G+T prior to surgical closure. The re-bleeding was controlled using a clip. This was deemed possibly related to the use of G+T.
- Subject had a re-bleeding of the target bleeding site treated with G+T prior to surgical closure. The re-bleeding was controlled using cautery. This was deemed possibly related to G+T and definitely related to the investigational procedure.

In sum, there was a single event identified as a SADE (subject who experienced skin reaction, described above).

Table 15: Serious Adverse Events

Note: Mortality events only occurred in the G+T control group 3/83 patients or 3.6% and serious adverse events tended to be related to re-bleeding episodes.

Adverse Event Type*	All	HEMOBLAST™	G+T
Abnormal Bloodwork	18/258 (7.0%)	13/175 (7.4%)	5/83 (6.0%)
Acute Kidney Injury	6/258 (2.3%)	4/175 (2.3%)	2/83 (2.4%)

Population	All	Site 2	Site 6	Site 8
EXPLORATORY ANALYSIS POPULATION*: Meet preoperative eligibility criteria but fail intraoperative eligibility criteria due to ineligible SBSS scores PLUS EFFICACY ANALYSIS POPULATION (excludes roll-in subjects)	25	9	7	8

* Four subjects were preoperatively eligible but then did not meet intraoperative criteria 02-04, 02-07, 02-08 and 08-03. Only subject 02-07 was excluded intraoperatively due to an SBSS score. Per the Statistical Analysis Plan, the exploratory analysis population should include the subjects intraoperatively ineligible resulting from one or both investigators providing a baseline SBSS score for the TBS that is in the ineligible range (i.e., SBSS 0, 4, or 5). Therefore, subject 02-07 was the only preoperatively eligible subject included in the exploratory analysis population.

Table 3: Number of Subjects in Each Analysis Population

Population	All	Abdominal	Orthopedic
TOTAL: Preoperative eligible + roll-in subjects	31	12	19
Intraoperative ineligible subjects	4	3	1
SAFETY ANALYSIS POPULATION: Preoperative AND Intraoperative eligible (ENROLLED and includes roll-in subjects)	27	9	18
Roll-in subjects	3	1	2
EFFICACY ANALYSIS POPULATION: Preoperative AND Intraoperative eligible (ENROLLED and EXCLUDES roll-in subjects)	24	8	16
EXPLORATORY ANALYSIS POPULATION: Meet pre-operative eligibility criteria but fail intra-operative eligibility criteria PLUS EFFICACY ANALYSIS POPULATION	25	9	16

C. Study Population Demographics and Baseline Parameters

The average age for all subjects was 62.8 and gender was roughly evenly distributed with 51.9% of all subjects being male and 48.1% female. In terms of race, the majority of subjects self-identified as Caucasian (77.8%) and the remaining subjects self-identified as African-American (22.2%). Self-identified Hispanic subjects were 22.2% of the entire subject pool. There were not any notable differences between sites in terms of demographic characteristics other than Site 2 enrolling all of the Hispanic subjects in the study while Sites 6 and 8 enrolled all of the African-American subjects in the study.

Table 4: Study Population Demographics

Measure	All (n=27)	Abdominal (n=9)	Orthopedic (n=18)
Age	62.8 ± 8.64 63.0 [55.0, 68.0]	60.4 ± 10.70 55.0 [53.0, 66.0]	64.0 ± 7.48 64.5 [62.0, 68.0]

Measure	All (n=27)	Abdominal (n=9)	Orthopedic (n=18)
Gender			
<i>Male</i>	14/27 (51.9%)	6/9 (66.7%)	8/18 (44.4%)
<i>Female</i>	13/27 (48.1%)	3/9 (33.3%)	10/18 (55.6%)
Ethnicity			
<i>Hispanic or Latino</i>	6/27 (22.2%)	6/9 (66.7%)	0/18 (0.0%)
<i>Not Hispanic or Latino</i>	21/27 (77.8%)	3/9 (33.3%)	18/18 (100.0%)
Race			
<i>Caucasian</i>	21/27 (77.8%)	9/9 (100.0%)	12/18 (66.7%)
<i>African American</i>	6/27 (22.2%)	0/9 (0.0%)	6/18 (33.3%)
<i>American Indian or Alaska Native</i>	0/27 (0.0%)	0/9 (0.0%)	0/18 (0.0%)
<i>Asian</i>	0/27 (0.0%)	0/9 (0.0%)	0/18 (0.0%)
<i>Native Hawaiian, Pacific Islander, or other</i>	0/27 (0.0%)	0/9 (0.0%)	0/18 (0.0%)
<i>Other</i>	0/27 (0.0%)	0/9 (0.0%)	0/18 (0.0%)

Table 5: Study Population Baseline Parameters

Measure	All	Abdominal	Orthopedic
Liver cell carcinoma	1/27 (3.7%)	1/9 (11.1%)	0/18 (0.0%)
Metastatic colon cancer to the liver	1/27 (3.7%)	1/9 (11.1%)	0/18 (0.0%)
Metastatic gastroesophageal junction cancer	1/27 (3.7%)	1/9 (11.1%)	0/18 (0.0%)
Metastatic rectal cancer	1/27 (3.7%)	1/9 (11.1%)	0/18 (0.0%)
Osteoarthritis of knee	10/27 (37.0%)	0/9 (0.0%)	10/18 (55.6%)
Pancreatic cancer	2/27 (7.4%)	2/9 (22.2%)	0/18 (0.0%)
Concomitant illnesses	20/27 (74.1%)	4/9 (44.4%)	16/18 (88.9%)
Surgical history related to the surgical area	9/27 (33.3%)	4/9 (44.4%)	5/18 (27.8%)
Diabetes			
<i>Type I</i>	0/27 (0.0%)	0/9 (0.0%)	0/18 (0.0%)
<i>Type II</i>	7/27 (25.9%)	4/9 (44.4%)	3/18 (16.7%)

Table 13 and Table 14 present the baseline SBSS score for the TBS by treatment group and surgical arm, respectively.

Table 13: Baseline SBSS Score for Each Treatment Group

Measure	All (N=242)	HEMOBLAST™ (N=159)	G+T (N=83)
Investigator SBSS Score			
0	0 (0.0%)	0 (0.0%)	0 (0.0%)
1	84 (34.7%)	59 (37.1%)	25 (30.1%)
2	100 (41.3%)	61 (38.4%)	39 (47.0%)
3	58 (24.0%)	39 (24.5%)	19 (22.9%)
4	0 (0.0%)	0 (0.0%)	0 (0.0%)
5	0 (0.0%)	0 (0.0%)	0 (0.0%)
SBSS score ≥2	158 (65.3%)	100 (62.9%)	58 (69.9%)

Numbers are n/N (percent).

Table 14: Baseline SBSS Score for Each Surgical Arm

Measure	All (N=242)	Cardiothoracic (N=58)	Abdominal (N=90)	Orthopedic (N=94)
Investigator SBSS Score				
0	0 (0.0%)	0 (0.0%)	0 (0.0%)	0 (0.0%)
1	84 (34.7%)	14 (24.1%)	8 (8.9%)	62 (66.0%)
2	100 (41.3%)	24 (41.4%)	46 (51.1%)	30 (31.9%)
3	58 (24.0%)	20 (34.5%)	36 (40.0%)	2 (2.1%)
4	0 (0.0%)	0 (0.0%)	0 (0.0%)	0 (0.0%)
5	0 (0.0%)	0 (0.0%)	0 (0.0%)	0 (0.0%)
SBSS score ≥2	158/242 (65.3%)	44/58 (75.9%)	82/90 (91.1%)	32/94 (34.0%)

Numbers are n/N (percent).

The orthopedic procedures had a much greater number of patients in the lower SBSS scores. The orthopedic procedures represented a great variety of procedures the largest numbers being in hip and knee replacement and other open surgeries. The abdominal procedure group of patients consisted of 90 patients. The most prevalent abdominal procedures included liver resection surgery for both primary and metastatic cancers with less than 5% of the liver resections performed for biliary cancers and benign hepatic tumors. There was a total of 32 liver resection cases of the 90 patients in the abdominal group. A large portion of the remaining abdominal surgeries included abdominoplasties and abdominal hernia repairs. The cardiothoracic cases comprised the smallest group of patients and consisted of patients requiring application of the subject or control hemostat to the sternal edges, aortotomy sites, saphenous vein graft anastomosis site (proximal or distal not defined), pericardium, and myocardial suture lines. Application to synthetic graft bleeding was excluded from this arm of the study.

Population	All	Cardiothoracic	Abdominal	Orthopedic
Lead-in subjects	16	6	7	3
EFFICACY ANALYSIS POPULATION: Preoperative AND Intraoperative eligible (ENROLLED and EXCLUDES lead-in subjects)	242 [∴]	58	90	94

* One subject failed to meet intra-operative eligibility criteria, due to an intraoperative SBSS score of 0, and one subject was withdrawn prior to hemostat application due to an intraoperative SBSS score of 0.

∴ Two subjects were enrolled, that were identified at subsequent monitoring visits, to have not met preoperative eligibility criteria. Because these subjects were enrolled and received HEMOBLAST™ or G+T application, they are considered part of the total, safety analysis, and efficacy analysis populations.

C. Study Population Demographics and Baseline Parameters

Table 12 presents physical measurement for height, weight, BMI, and blood pressure, which were similar between treatment arms.

Table 12: Physical Measurements by Treatment Group

Measure	All (n=258)	HEMOBLAST™ (n=175)	G+T (n=83)
Height (cm)	167.4 ± 10.19	167.4 ± 10.43	167.3 ± 9.72
	167.0 [160.0, 175.0]	165.0 [160.0, 175.0]	167.0 [160.0, 175.0]
Weight (kg)	83.8 ± 19.121	84.6 ± 19.19	82.1 ± 19.26
	82.5 [70.5, 97.0]	83.0 [71.0, 97.5]	80.9 [69.4, 97.0]
BMI (kg/m ²)	29.9 ± 6.33	30.2 ± 6.14	29.3 ± 6.71
	29.1 [25.5, 33.7]	29.4 [25.6, 34.1]	27.9 [25.1, 33.1]
Systolic Blood pressure (mmHg)	129.0 ± 18.54	128.8 ± 18.28	129.4 ± 19.18
	128.0 [116.0, 138.0]	128.0 [116.0, 138.0]	128.0 [118.0, 141.0]
Diastolic Blood pressure (mmHg)	76.6 ± 12.60	76.5 ± 12.33	76.9 ± 13.24
	78.0 [70.0, 85.0]	78.0 [69.0, 84.0]	76.0 [70.0, 85.0]

Numbers are mean ± SD (N)/ median.

Measure	All	Abdominal	Orthopedic
Smoking history			
<i>Never</i>	15/27 (55.6%)	4/9 (44.4%)	11/18 (61.1%)
<i>Former smoker</i>	9/27 (33.3%)	4/9 (44.4%)	5/18 (27.8%)
<i>Occasionally or socially</i>	0/27 (0.0%)	0/9 (0.0%)	0/18 (0.0%)
<i>Less than half a pack a day</i>	0/27 (0.0%)	0/9 (0.0%)	0/18 (0.0%)
<i>Between half a pack and one pack a day</i>	3/27 (11.1%)	1/9 (11.1%)	2/18 (11.1%)
<i>More than one pack a day</i>	0/27 (0.0%)	0/9 (0.0%)	0/18 (0.0%)
Malignancies and prior therapies	10/27 (37.0%)	6/9 (66.7%)	4/18 (22.2%)
Hepatic Disease	3/27 (11.1%)	3/9 (33.3%)	0/18 (0.0%)

Numbers are n/N (percent) for categorical measures.

Table 6: Study Population Surgical Baseline Parameters

Measure	All	Site 2	Site 6	Site 8	Abdominal	Orthopedic
Surgical Indication						
<i>Abdominal</i>	9/27 (33.3%)	9/9 (100.0%)	0/8 (0.0%)	0/10 (0.0%)	9/9 (100.0%)	0/18 (0.0%)
<i>Orthopedic</i>	18/27 (66.7%)	0/9 (0.0%)	8/8 (100.0%)	10/10 (100.0%)	0/9 (0.0%)	18/18 (100.0%)
Location						
<i>Medial Retinaculum</i>	1/27 (3.7%)	0/9 (0.0%)	0/8 (0.0%)	1/10 (10.0%)	0/9 (0.0%)	1/18 (5.6%)
<i>Quad</i>	3/27 (11.1%)	0/9 (0.0%)	0/8 (0.0%)	3/10 (30.0%)	0/9 (0.0%)	3/18 (16.7%)
<i>Quad Tendon</i>	3/27 (11.1%)	0/9 (0.0%)	0/8 (0.0%)	3/10 (30.0%)	0/9 (0.0%)	3/18 (16.7%)
<i>Segment 2</i>	1/27 (3.7%)	1/9 (11.1%)	0/8 (0.0%)	0/10 (0.0%)	1/9 (11.1%)	0/18 (0.0%)
<i>Segment 3</i>	3/27 (11.1%)	3/9 (33.3%)	0/8 (0.0%)	0/10 (0.0%)	3/9 (33.3%)	0/18 (0.0%)
<i>Segment 4</i>	2/27 (7.4%)	2/9 (22.2%)	0/8 (0.0%)	0/10 (0.0%)	2/9 (22.2%)	0/18 (0.0%)
<i>Segment 5</i>	2/27 (7.4%)	2/9 (22.2%)	0/8 (0.0%)	0/10 (0.0%)	2/9 (22.2%)	0/18 (0.0%)
<i>Segments 4 & 6</i>	1/27 (3.7%)	1/9 (11.1%)	0/8 (0.0%)	0/10 (0.0%)	1/9 (11.1%)	0/18 (0.0%)
<i>Suprapatellar pouch</i>	10/27 (37.0%)	0/9 (0.0%)	8/8 (100.0%)	2/10 (20.0%)	0/9 (0.0%)	10/18 (55.6%)
<i>Tendon</i>	1/27 (3.7%)	0/9 (0.0%)	0/8 (0.0%)	1/10 (10.0%)	0/9 (0.0%)	1/18 (5.6%)

Measure	All	Site 2	Site 6	Site 8	Abdominal	Orthopedic
Tissue type						
Bone	1/27 (3.7%)	0/9 (0.0%)	1/8 (12.5%)	0/10 (0.0%)	0/9 (0.0%)	1/18 (5.6%)
Liver	9/27 (33.3%)	9/9 (100.0%)	0/8 (0.0%)	0/10 (0.0%)	9/9 (100.0%)	0/18 (0.0%)
Muscle And Soft Tissue	3/27 (11.1%)	0/9 (0.0%)	3/8 (37.5%)	0/10 (0.0%)	0/9 (0.0%)	3/18 (16.7%)
Quad/Tendon	1/27 (3.7%)	0/9 (0.0%)	0/8 (0.0%)	1/10 (10.0%)	0/9 (0.0%)	1/18 (5.6%)
Subcutaneous soft tissue	3/27 (11.1%)	0/9 (0.0%)	0/8 (0.0%)	3/10 (30.0%)	0/9 (0.0%)	3/18 (16.7%)
Synovium	2/27 (7.4%)	0/9 (0.0%)	2/8 (25.0%)	0/10 (0.0%)	0/9 (0.0%)	2/18 (11.1%)
Synovium And Muscle Tissue	2/27 (7.4%)	0/9 (0.0%)	2/8 (25.0%)	0/10 (0.0%)	0/9 (0.0%)	2/18 (11.1%)
Tendon	5/27 (18.5%)	0/9 (0.0%)	0/8 (0.0%)	5/10 (50.0%)	0/9 (0.0%)	5/18 (27.8%)
Tendon/Soft Tissue	1/27 (3.7%)	0/9 (0.0%)	0/8 (0.0%)	1/10 (10.0%)	0/9 (0.0%)	1/18 (5.6%)
TBS approximate dimensions (cm²)	21.1 ± 75.88 (27) 4.0 [4.0, 9.0]	51.0 ± 131.04 (9) 4.0 [4.0, 16.0]	3.7 ± 1.53 (8) 4.0 [3.1, 4.0]	8.1 ± 4.07 (10) 8.3 [4.0, 12.0]	51.0 ± 131.04 (9) 4.0 [4.0, 16.0]	6.1 ± 3.84 (18) 4.0 [4.0, 9.0]
Conventional Procedures for Hemostasis						
Pressure	12/27 (44.4%)	0/9 (0.0%)	2/8 (25.0%)	10/10 (100.0%)	0/9 (0.0%)	12/18 (66.7%)
Suture, ligation	0/27 (0.0%)	0/9 (0.0%)	0/8 (0.0%)	0/10 (0.0%)	0/9 (0.0%)	0/18 (0.0%)
Cautery	2/27 (7.4%)	0/9 (0.0%)	1/8 (12.5%)	1/10 (10.0%)	0/9 (0.0%)	2/18 (11.1%)
Other	0/27 (0.0%)	0/9 (0.0%)	0/8 (0.0%)	0/10 (0.0%)	0/9 (0.0%)	0/18 (0.0%)
None (impractical)	14/27 (51.9%)	9/9 (100.0%)	5/8 (62.5%)	0/10 (0.0%)	9/9 (100.0%)	5/18 (27.8%)

Numbers are mean ± SD (N)/ median for continuous measures, and n/N (percent) for categorical measures.

D. Safety and Effectiveness Results

The mean paired kappa across all sites was 0.9301 indicative of almost perfect agreement. Given the pre-specified criteria to advance to the pivotal clinical investigation was that the mean paired kappa statistic for assignment of SBSS scores was >0.80, this endpoint was deemed successful. These results validated the use of the SBSS for the pivotal study.

There were 41 minor adverse events not believed to be device related. One serious adverse event, portal vein thrombosis, occurred in a patient who underwent a right hepatic lobectomy for metastatic colon cancer with metastectomy of additional metastases in the left hepatic lobe. This required prolonged operative time and Pringle maneuver all of which more likely contributed to the portal vein thrombosis. This event was not likely caused by embolization of Hemoblast applied to the liver parenchymal

B. Accountability of PMA Cohort

The study was stopped early for efficacy, per the IDMC recommendation, based on the pre-specified stopping rules. At the time of the interim analysis, a total of 258 subjects were enrolled in the study, including 16 lead-in subjects and 242 randomized subjects. The interim analysis on completion of half the intended total number of patients in the pivotal study demonstrated non-inferiority at the 6 minute hemostasis primary endpoint and superiority at the 10-30% superiority at all the secondary endpoints including hemostasis at 6 minutes, 3 minutes, and 10 minutes as well as superiority on time to preparation for use. Furthermore, the pivotal study demonstrated non inferiority of Hemoblast for hemostasis at 10 minutes.

Table 10 represents the enrollment for each treatment group and **Table 11** represents the enrollment for each surgical arm. Randomization was stratified by surgical arm.

Table 10: Enrollment Details for Each Treatment Group

Population	All	HEMOBLAST™	G+T
TOTAL: Preoperative eligible + lead-in subjects	260* [‡]	175	83
Intraoperative ineligible subjects	1	0	0
SAFETY ANALYSIS POPULATION: Preoperative AND Intraoperative eligible (ENROLLED and includes lead-in subjects)	258 [‡]	175	83
Lead-in subjects	16	16	0
EFFICACY ANALYSIS POPULATION: Preoperative AND Intraoperative eligible (ENROLLED and EXCLUDES lead-in subjects)	242 [‡]	159	83

* One subject failed to meet intra-operative eligibility criteria, due to an intraoperative SBSS score of 0, and one subject was withdrawn after randomization/enrollment, but prior to hemostat application.

[‡] Two subjects were enrolled, that were identified at subsequent monitoring visits, to have not met preoperative eligibility criteria. Because these subjects were enrolled and received HEMOBLAST™ or G+T application, they are considered part of the total, safety analysis, and efficacy analysis populations.

Table 11: Enrollment by Surgical Arm

Population	All	Cardiothoracic	Abdominal	Orthopedic
TOTAL: Preoperative eligible + lead-in subjects	260* [‡]	64	98	98
Intraoperative ineligible subjects	1	0	0	1
SAFETY ANALYSIS POPULATION: Preoperative AND Intraoperative eligible (ENROLLED and includes lead-in subjects)	258 [‡]	64	97	97

2. Patient Follow-up Schedule

All patients underwent the same intraoperative investigational evaluations. During the surgery, hemostasis was evaluated by the investigator 3, 6, and 10 minutes after application of the investigational or control treatment until hemostasis was achieved. Reapplication of the randomized hemostat was performed at the 3 and 6 minute evaluation time points, as needed. In cases where hemostasis was not achieved by 10 minutes, the Investigator may have used whatever means necessary in order to control bleeding, except for any hemostatic products containing thrombin or aprotinin. Thrombin should not have been used in subjects randomized to receive HEMOBLAST, but may have been used in subjects randomized to the G+T arm. Investigators were permitted to use any of the remaining randomized hemostat for bleeding sites other than the target bleeding site (“TBS”) or any hemostatic product not containing thrombin or aprotinin.

For cardiothoracic procedures, the TBS was identified after heparin reversal. Additionally, the TBS was not at the site of a valve implant or at the site of a synthetic graft or patch implant. Bleeding severity and hemostasis were assessed using the SBSS.

Safety assessments occurred one day and 6 weeks postoperatively. Blood draws for antibody evaluation were performed preoperatively (within 4 weeks of surgery) and 6 weeks postoperatively.

3. Clinical Endpoints

The primary efficacy endpoint of this clinical investigation was non-inferiority of HEMOBLAST™ relative to G+T for success at achieving hemostasis within 6 minutes.

The secondary efficacy endpoints of this clinical investigation were:

- Superiority of HEMOBLAST™ relative to G+T in mean preparation time from the opening of package to product being ready to use;
- Non-inferiority of HEMOBLAST™ relative to G+T for success at achieving hemostasis within 3 minutes;
- Superiority of HEMOBLAST™ relative to G+T for success at achieving hemostasis within 6 minutes; and
- Superiority of HEMOBLAST™ relative to G+T for success at achieving hemostasis within 3 minutes.

The pivotal study was designed to test the non-inferiority hypothesis of HEMOBLAST™ relative to G+T for success at achieving hemostasis within 6 minutes and within a 10% margin. The study pre-specified a single interim analysis for early stopping due to futility or effectiveness after outcome data had been observed on 240 treated subjects. If the study were to continue to the final analysis with a decision in favor of comparable effectiveness, it was anticipated that efficacy data would be available on a maximum of 400 treated subjects (approximately 267 patients treated with HEMOBLAST™ under a 2:1 randomization scheme).

Effectiveness analyses were conducted on the time-to-hemostasis (“TTH”) population, defined as all subjects who were randomized, received study intervention, and had a TTH assessment recorded regardless of whether the measurement was censored (defined as the use of an additional hemostatic product or surgical rescue prior to the end of observation time, or failing to achieve and maintain complete hemostasis prior to the end of observation time). Lead-in subjects were not part of the TTH Population. Safety analyses were conducted on the Full Analysis population, defined as all subjects who were randomized into the study and received study intervention and all lead-in subjects.

In the primary analyses, missing TTH values were not imputed. All values right censored prior to 6 minutes were considered treatment failures for the purpose of the primary analysis.

edge of resection and resolved with anticoagulation.

Table 7: Adverse Events

Measure	All	Site 2	Site 6	Site 8	Abdominal	Orthopedic
Number of AEs	41	29	3	9	29	12
Number of subjects with AEs	15/27 55.6% (37.3%,72.4%)	8/9 88.9% (56.5%, 98.0%)	2/8 25.0% (7.1%, 59.1%)	5/10 50.0% (23.7%, 76.3%)	8/9 88.9% (56.5%, 98.0%)	7/18 38.9% (20.3%, 61.4%)
Number of AEs by severity						
<i>Mild</i>	16	13	0	3	13	3
<i>Moderate</i>	22	13	3	6	13	9
<i>Severe</i>	3	3	0	0	3	0
Number of subjects with AEs by severity*						
<i>Mild</i>	5/27 18.5% (8.2%, 36.7%)	4/9 44.4% (18.9%, 73.3%)	0/8 0.0% (0.0%, 32.4%)	1/10 10.0% (1.8%, 40.4%)	4/9 44.4% (18.9%, 73.3%)	1/18 5.6% (1.0%, 25.8%)
<i>Moderate</i>	8/27 29.6% (15.9%, 48.5%)	2/9 22.2% (6.3%, 54.7%)	2/8 25.0% (7.1%, 59.1%)	4/10 40.0% (16.8%, 68.7%)	2/9 22.2% (6.3%, 54.7%)	6/18 33.3% (16.3%, 56.3%)
<i>Severe</i>	2/27 7.4% (2.1%, 23.4%)	2/9 22.2% (6.3%, 54.7%)	0/8 0.0% (0.0%, 32.4%)	0/10 0.0% (0.0%, 27.8%)	2/9 22.2% (6.3%, 54.7%)	0/18 0.0% (0.0%, 17.6%)
Number AEs related to device	0	0	0	0	0	0
Number of subjects with AEs related to device	0/27 0.0% (0.0%, 12.5%)	0/9 0.0% (0.0%, 29.9%)	0/8 0.0% (0.0%, 32.4%)	0/10 0.0% (0.0%, 27.8%)	0/9 0.0% (0.0%, 29.9%)	0/18 0.0% (0.0%, 17.6%)
Number of SAEs	8	7	1	0	7	1
Number of subjects with SAEs	4/27 14.8% (5.9%, 32.5%)	3/9 33.3% (12.1%, 64.6%)	1/8 12.5% (2.2%, 47.1%)	0/10 0.0% (0.0%, 27.8%)	3/9 33.3% (12.1%, 64.6%)	1/18 5.6% (1.0%, 25.8%)
Number of SADEs	0	0	0	0	0	0
Number of subjects with SADEs	0/27 0.0% (0.0%, 12.5%)	0/9 0.0% (0.0%, 29.9%)	0/8 0.0% (0.0%, 32.4%)	0/10 0.0% (0.0%, 27.8%)	0/9 0.0% (0.0%, 29.9%)	0/18 0.0% (0.0%, 17.6%)
Number of UADEs	0	0	0	0	0	0
Number of subjects with UADEs	0/27 0.0% (0.0%, 12.5%)	0/9 0.0% (0.0%, 29.9%)	0/8 0.0% (0.0%, 32.4%)	0/10 0.0% (0.0%, 27.8%)	0/9 0.0% (0.0%, 29.9%)	0/18 0.0% (0.0%, 17.6%)

Numbers are n/N percent (95% CI).

Wilson confidence limits (score based) are used in the table.

*The most severe AE is counted for each subject.

Table 8: Thromboembolic Events

	Total Thromboembolic Events	Pulmonary Embolism	Deep Venous Thrombosis	Stroke	Thrombosed AV Fistula	Other(portal vein thrombosis)
HEMOBLAST™ US Pilot	1/27 (3.7%)	0/27 (0.0%)	0/27 (0.0%)	0/27 (0.0%)	0/27 (0.0%)	1/27 (3.7%)

HEMOBLAST™ induced hemostasis at 3 minutes for half of all subjects in the efficacy analysis population and 79.2% of this group (19/24) achieved hemostasis at 6 minutes. There were 2/24 subjects that did not achieve hemostasis at 10 minutes. Comparing abdominal surgery subjects to orthopedic surgery subjects at the 6-minute time point, it appears the orthopedic group performed better with 93.8% of subjects achieving hemostasis compared to 50.0% of subjects in the abdominal surgery group.

Table 9: Proportion of Subjects Achieving Hemostasis at 3, 6, and 10 Minutes - Efficacy Analysis Population

Time	All	Site 2	Site 6	Site 8	Abdominal	Orthopedic
3 minutes	12/24 50.0% (31.4%, 68.6%)	2/8 25.0% (7.1%, 59.1%)	3/7 42.9% (15.8%, 75.0%)	7/9 77.8% (45.3%, 93.7%)	2/8 25.0% (7.1%, 59.1%)	10/16 62.5% (38.6%, 81.5%)
6 minutes	19/24 79.2% (59.5%, 90.8%)	4/8 50.0% (21.5%, 78.5%)	6/7 85.7% (48.7%, 97.4%)	9/9 100.0% (70.1%, 100.0%)	4/8 50.0% (21.5%, 78.5%)	15/16 93.8% (71.7%, 98.9%)
10 minutes	22/24 91.7% (74.2%, 97.7%)	6/8 75.0% (40.9%, 92.9%)	7/7 100.0% (64.6%, 100.0%)	9/9 100.0% (70.1%, 100.0%)	6/8 75.0% (40.9%, 92.9%)	16/16 100.0% (80.6%, 100.0%)

Numbers are n/N percent (95% CI).

Wilson confidence limits (score based) are used in the table.

Cumulative numbers of subjects achieving hemostasis at each time are counted.

Pivotal Study:

A. Study design

This was a prospective, randomized, controlled, multicenter pivotal study to evaluate the safety and efficacy of the HEMOBLAST™ Bellows (“HEMOBLAST™”) in cardiothoracic, abdominal (both soft tissue and organ space), and orthopedic lower extremity surgeries. Subjects were randomized intraoperatively in a 2:1 ratio to HEMOBLAST™ or the control treatment, an absorbable gelatin sponge, USP and recombinant thrombin (“G+T”). The study also included a lead-in phase, to ensure correct device application procedures, wherein subjects were not randomized. Lead-in subjects received HEMOBLAST™ and were followed for safety only.

The specific objective of the study was to evaluate the safety and efficacy of HEMOBLAST™ compared to G+T. The primary hypothesis in this study was that HEMOBLAST™ is non-inferior relative to G+T for success at achieving hemostasis within 6 minutes. A Clinical Events Committee (“CEC”) and Independent Data Monitoring Committee (“IDMC”) were utilized to review all serious adverse events and safety and efficacy data, respectively.

1. Clinical Inclusion and Exclusion Criteria

A subject must have met all the following preoperative inclusion criteria to be enrolled into the investigation:

- Subject is undergoing an elective open cardiothoracic, abdominal, or orthopedic lower extremity surgery;
- Subject or an authorized legal representative is willing and able to give prior written informed consents for investigation participation;
- Subject undergoing cardiothoracic surgery is not allergic to protamine; and
- Subject is 21 years of age or older.

A subject must have met all the following intraoperative inclusion criteria to be enrolled into the investigation:

- Subject does not have an active or suspected infection at the surgical site;
- Subject undergoing cardiothoracic surgery with anticoagulation must have anticoagulation reversed prior to Target Bleeding Site (TBS) identification and treatment;
- Subject in whom the Investigator is able to identify a TBS for which any applicable conventional means for hemostasis are ineffective or impractical; and
- Subject has a TBS with an SBSS score of 1, 2, or 3.

Subjects were excluded if they met any of the following exclusion criteria:

- Subject is undergoing a laparoscopic, thoracoscopic, or robotic surgical procedure;
- Subject is undergoing a neurologic surgical procedure;
- Subject is undergoing a spinal surgical procedure;
- Subject is undergoing an emergency surgical procedure;
- Subject is pregnant, planning on becoming pregnant during the follow-up period, or actively breast-feeding;
- Subject has a clinically significant coagulation disorder or disease, defined as a platelet count < 100,000 per microliter or International Normalized Ratio > 1.5 within 4 weeks of surgery;
- Subject receiving intravenous heparin within 12 hours before surgery or oral Coumadin within 2 days before surgery;
- Subject receiving antiplatelet medications within 5 days prior to surgery;
- Subject undergoing abdominal or orthopedic lower extremity surgery receiving aspirin within 7 days prior to surgery;
- Subject has an active or suspected infection at the surgical site;
- Subject has had or has planned to receive any organ transplantation;
- Subject has a known sensitivity or allergy to bovine and/or porcine substance(s) or any other component(s) of the hemostatic agent;
- Subject has ASA classification of 5;
- Subject has a life expectancy of less than 3 months;
- Subject has a known psychiatric disorder, which in the opinion of the Principal Investigator, would preclude the subject from completing this clinical study;
- Subject has a documented severe congenital or acquired immunodeficiency;
- Subject has religious or other objections to porcine, bovine, or human components;
- Subject in whom the investigational or control device will be used at the site of a valve replacement or repair;
- Subject in whom the investigational or control device will be used at the site of a synthetic graft or patch implant;
- Subject is currently participating or has participated in another clinical trial within the past 30 days and is receiving/has received an investigational drug, device, or biologic agent; and
- Subject is not appropriate for inclusion in the clinical trial, per the medical opinion of the Principal Investigator.

HEMOBLAST™ Bellows

Laparoscopic Applicator

Caution: Federal law restricts this device to sale on or by the order of a physician.














This device is intended for use with HEMOBLAST™ Bellows. Refer to HEMOBLAST™ Bellows Instructions for Use for complete information.

1. Device Description	p.3	5. Directions for Use	p.3
2. How Supplied	p.3	6. Storage and Handling	p.4
3. Intended Use / Indications for Use	p.3	7. Sterilization	p.4
4. Warnings and precautions	p.3	8. Symbols	p.6



Manufacturer
Dilon Technologies France SAS
8, Allée Irène Joliot-Curie St-Priest
69800, France
www.dilon.com
Phone: +1 877 GO DILON
Fax: +1 757 269 4912

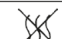8. Symbols

Symbol	Standard	Reference #	Title	Definition
	ISO 15223-1	5.1.5	Batch code	Indicates the manufacturer's batch code so that the batch or lot can be identified.
	ISO 15223-1	5.1.6	Catalogue Number	Indicates the manufacturer's catalogue number so that the medical device can be identified.
	ISO 15223-1	5.1.4	Use-By date	Indicates the date after which the medical device is not to be used.
	ISO 15223-1	5.3.7	Temperature Limit Store between 2°C and 25°C (36-77°F)	Indicates the temperature limits to which the medical device can be safely exposed.
	ISO 15223-1	5.2.6	Do not resterilize	Indicates a medical device that is not to be resterilized.
	ISO 15223-1	5.2.8	Do not use if package is damaged and consult instructions for use	Indicates a medical device that should not be used if the package has been damaged or opened and that the user should consult the instructions for use for additional information.
	ISO 15223-1	5.4.2	Do not re-use	Indicates a medical device that is intended for one single use only.
	ISO 15223-1	5.4.3	Consult the instruction for use	Indicates the need for the user to consult the instructions for use.
	ISO 15223-1	5.2.4	Sterilized using irradiation	Indicates a medical device that has been sterilized using irradiation.
	ISO 15223-1	5.2.12	Double sterile barrier system	Indicates a double sterile barrier system.
	ISO 15223-1	5.6.3	Non-pyrogenic	Indicates a medical device that is non-pyrogenic.
	FDA Guidance "Alternative to Certain Prescription Device Labeling Requirements", issued 1/21/2000	N/A	Caution: Federal law restricts this device to sale on or by the order of a physician.	Requires prescription in the United States.
	ISO 15223-1	5.1.1	Manufacturer	Indicates the medical device Manufacturer

IFU-LAP01-US-V02 2023-02-22

1. Device Description

The HEMOBLAST™ Bellows Laparoscopic Applicator is a sterile, single use laparoscopic applicator device. The HEMOBLAST™ Bellows Laparoscopic Applicator is designed to be used with HEMOBLAST™ Bellows to deliver the hemostatic powder during laparoscopic surgical procedures. The laparoscopic applicator is 35 cm in length. HEMOBLAST™ Bellows is available separately.

2. How Supplied

The HEMOBLAST™ Bellows Laparoscopic Applicator is supplied sterile in double packaging. One (1) box contains six (6) individually packaged laparoscopic applicators.

3. Indications

The HEMOBLAST™ Bellows Laparoscopic Applicator is indicated for use in delivering HEMOBLAST™ Bellows hemostatic powder to bleeding surgical sites through a 5 mm or larger diameter trocar.

4. Warnings and Precautions

- Do not use the laparoscopic applicator to manipulate or retract organs or tissue.
- To prevent product contamination prior to application, always follow aseptic techniques.
- The product must be manipulated and used by qualified personnel according to the general principles of sterility and pre-medication.
- The tip of the laparoscopic applicator must be directly visualized at all times to minimize potential for unintended contact with tissue, organs, or blood as well as any possible unintended spillage of residual powder.
- The laparoscopic applicator should not be bent over (>) 65° to avoid any kinking risk.
- The safety and effectiveness of the combined use of the laparoscopic applicator with other hemostatic agents has not been evaluated and is therefore not recommended.
- Ensure all gauze and other accessories placed in the abdomen are removed after use.

5. Directions for Use

Applicator Assembly

1. Inspect the integrity of the product and its packaging prior to use. Do not use if damaged.
2. Peel open the bags to remove the laparoscopic applicator from its packaging with attention to sterile procedures.
3. Ensure that both ends of the laparoscopic applicator are rounded and not sharp.
4. Remove the HEMOBLAST™ Bellows from its packaging and remove its cap using a twisting motion. Discard cap.
5. Firmly insert the laparoscopic applicator into the nozzle of HEMOBLAST™ Bellows until resistance can be felt and a firm connection is observed. Either end of the laparoscopic applicator can be inserted into the HEMOBLAST™ Bellows nozzle.
6. The assembled device can then be used for a laparoscopic procedure.

Clinical Application

1. Prepare the application site as outlined in the HEMOBLAST™ Bellows Instructions for Use (IFU).
2. Place one or more dry gauze in close proximity of the target bleeding site using laparoscopic forceps.
3. Place one or more saline-soaked gauze(s) in close proximity of the target bleeding site using laparoscopic forceps (avoid contact with blood) prior to applying HEMOBLAST™ Bellows hemostatic powder on the wound.
4. Gently insert the laparoscopic applicator through a 5 mm or larger diameter trocar.
5. Avoid applying excessive bending forces on the laparoscopic applicator to avoid kinking. If kinking occurs, do not use the laparoscopic applicator; remove it and replace with a new laparoscopic applicator.
6. Avoid direct contact between the laparoscopic applicator tip and blood at the surgical site to minimize the potential for clogging the tip during HEMOBLAST™ Bellows application.
7. If the laparoscopic applicator becomes occluded/clogged during use, replace with a new laparoscopic applicator. Do not cut the laparoscopic applicator to remove the occlusion.
8. Use dry gauze to blot excess blood from the target bleeding site so the hemostatic powder may be applied directly to the source of bleeding. The wound surface should be as dry as possible before application of HEMOBLAST™ Bellows hemostatic powder.
9. Immediately position the laparoscopic applicator approximately 2-3 cm from the site of bleeding. Laparoscopic graspers may be used to facilitate precise application.
10. Apply the HEMOBLAST™ Bellows powder as directed in the HEMOBLAST™ Bellows IFU.
11. When applied with the laparoscopic applicator tip 2-3 cm from the target bleeding site, HEMOBLAST™ Bellows can treat a surface area up to 50 cm².
12. Immediately use a saline-soaked gauze to hold the hemostatic powder in contact with the bleeding area as directed in the HEMOBLAST™ Bellows IFU.
13. Gently remove the laparoscopic applicator through the trocar.

6. Storage and Handling

- The product is sterile unless the package has been opened, damaged, or otherwise contaminated.
- For single-use only.
- Do not reuse or resterilize. Reuse of single-use devices creates a potential risk of patient or user infections.
- Do not use after the expiration date indicated on the package.
- Product must be stored in a dry environment at room temperature.
- The product must be discarded after use according to the institutional policy for medical devices.

7. Sterilization

HEMOBLAST™ Bellows Laparoscopic Applicator is sterilized by gamma irradiation.

For any additional information, please contact the manufacturer.